

 **HEALFASTER**

G A M E  R E A D Y[®]
ACTIVATE NEXT-LEVEL INJURY & POST-OP RECOVERY



Heal Faster with the Game Ready

Accelerated Recovery System

Heal Faster is dedicated to helping people get back to their best after injury or surgery. Whether you are a post-op surgical patient or an elite athlete with a sports injury, the Game Ready system helps you get back into action as quickly as possible.

What is Game Ready?

Game Ready is the ONLY device that uses a combination of cold and compression therapies to accelerate recovery following injury or surgery.

Game Ready uses a combination of anatomically engineered wraps attached to a control unit that is filled with ice and water to deliver therapeutic cryotherapy (cold) consistently to the injury or surgical site. At the same time, the Game Ready system progressively compresses and releases pressure to the targeted area.

The device mimics the body's natural muscle contractions while cooling the tissue, helping to proactively aid lymphatic function and encourage cellular oxygen. For you that means less pain, less medication and a faster return to the activities you love.

Game Ready is so effective, it is the number one recovery technology for patients, elite athletes, athletic trainers, professional teams in every sport, special military forces, and the world's foremost orthopaedic surgeons, sports medicine physicians, and physical therapists.

Patients love the Game Ready system because it is intuitive and easy to use.

Game Ready's control unit allows you to create custom settings with a specific temperature, pressure level, and treatment time. There are also preset programs available so you can easily and proactively control your recovery with the touch of a button.

“Patients love the Game Ready system because it is intuitive and easy to use.”



Game Ready is the **international leader** in innovative systems for sports injury recovery and post-op orthopaedic rehabilitation.

Human Inspired, Technology Driven

The combination of cold and compression therapies allows you to reduce swelling and pain without narcotics or pain medication.

Our cryotherapy wraps have integrated chambers that circulate cold water and pneumatic air in separate, dedicated loop systems. The cold water flows back through an ice reservoir in the control unit to maintain the desired cold temperature.



At the same time, the pneumatic air provides a pumping effect that increases the flow of fluids in the desired area. This action helps pump away excess edema to reduce swelling and increases blood flow to promote healing. Intermittent compression enhances the benefits of therapeutic cold by allowing it to penetrate deeper and last longer than other methods.

CRYOTHERAPY (COLD)

Reduces Swelling	Reduces inflammation and swelling
Reduces Tissue Damage	Reduces metabolic activity and cellular demand for oxygen, reducing secondary tissue damage
Temperature Reduction	Draws heat away from the injury site for more efficient temperature reduction


COMPRESSION

Edema Prevention	Aids in the prevention of edema formation and removal of swelling
Reduces Swelling	Mimics natural muscle contractions and “pumps” edema/swelling out of the injured area and reduces swelling
Optimizes Drainage	Optimizes lymphatic drainage
Increases Blood Flow	Increases blood flow and oxygen delivery
Faster Healing Process	Stimulates tissue repair and encourages acceleration of the healing process

Why choose Game Ready?



“ Compared to traditional RICE (Rest-Ice-Compression-Elevation) applications, the Game Ready System provides greater therapeutic benefits.

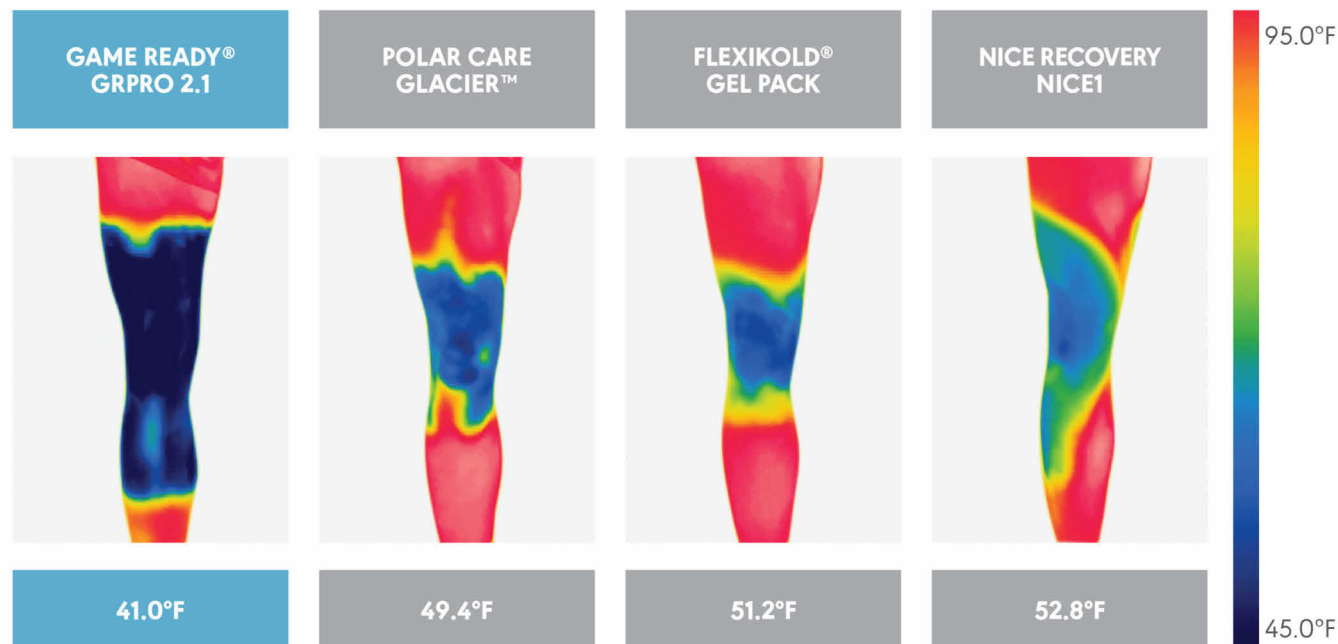

Thermal images and minimum skin temperatures taken immediately after 30 minutes of treatment.

The infrared thermal imagery above shows that Game Ready provides faster, colder, and longer-lasting therapy when compared to other products. The blue area indicates how Game Ready cools better and covers a larger area enabling you to heal faster.

As an added benefit, Game Ready's active compression enhances recovery when compared to using cold

therapy alone. This is because Game Ready's anatomic, circumferential wraps envelope the injury or surgery site and then conform to the contours of the body with pneumatic compression, assuring better surface contact and more effective cooling.

“ Game Ready has been proven to reduce narcotic consumption by 25%.



Over **9700 patients** were surveyed about Game Ready:

98%

were glad their physician offered Game Ready as an option for recovery.

87%

said Game Ready helped reduce the amount of pain medication they took, or enabled them to stop taking pain medication sooner.

96%

said Game Ready provided a better post-op recovery than a previous surgery using a different cold therapy.

Game Ready Wraps for Head-to-Toe Coverage

Game Ready wraps come in a variety of sizes and fits that are easy to apply and comfortable to wear. You can even wear them overnight while sleeping!

No matter what your injury or surgery, Game Ready has a wrap that can help accelerate your recovery.

Better **coverage**. Better **fit**. Better **results**.

With conventional cryotherapy such as ice and gel packs, the cold source lays statically over the injury site, quickly warming and becoming ineffective as it absorbs the body's heat. Game Ready wraps actively draw heat away from the injury by continuously replenishing cold water drawn from the ice reservoir.

The water that is continually circulating through the wrap is then returned back to the ice reservoir, cooled, then recirculated again. This ensures the continuous cooling effect needed to heal faster.



How to Order

Choose from one of the options below:

1.

Reserve your Game Ready through your **healthcare provider**. They will take payment and make the necessary arrangements.

2.

Reserve your Game Ready online at **www.healfaster.ca** with a simple to complete form.

3.

Call our toll free number **1-844-441-FAST (3278)** and one of our team will be happy to help you complete your order.

\$250.00 + Tax Per Week

2-Week Rental

Free Shipping Included

Game Ready FAQ

Our website contains complete how-to videos and user guides, and our 24/7 customer support service ensures that you know how to easily operate the Game Ready accelerated recovery system.

1. How much does Game Ready cost?

You can rent a Game Ready device for \$250 + tax per week with a 2-week minimum rental. This includes the cost of return shipping.

2. Is Game Ready covered by insurance?

Some insurance providers will cover the cost of the Game Ready rental. Please contact your provider for details. Provincial workers' compensation boards will cover the cost up to an 8-week rental.

3. How long does it take to order a Game Ready?

You can rent a Game Ready device in minutes over the phone, on our website, or through your healthcare provider. We can ship the unit to you the same day or next day. We strive to ensure you have your device 2 days prior to your surgery so that you can feel secure knowing you know how to use the Game Ready when you need it most.

4. How long do I need to use Game Ready for?

Your physician will tell you how long you need to use Game Ready. Once you have the device we believe "more is better" - the longer you use Game Ready each day the better you will feel and the quicker you will recover.

G A M E R E A D Y[®]

ACTIVATE NEXT-LEVEL INJURY & POST-OP RECOVERY

Heal Faster is proud to be the Canadian
Distributor of Game Ready products

 **HEALFASTER**

Heal Faster Canada
170 Zenway Blvd., Unit #3
Woodbridge, ON L4H 2Y7
1-844-441-FAST (3278)
info@healfaster.ca
www.healfaster.ca



We Offer 24/7 Customer Support Service